

Real-world clinical performance of commercial SARS-CoV-2 rapid antigen tests in suspected COVID-19: A systematic meta-analysis of available data as per November 20, 2020

Johannes Hayer¹, Dusanka Kasapic², Claudia Zemmrich^{2,3}

¹Roche Diagnostics GmbH, Mannheim, Germany

²Roche Diagnostics International Ltd. Rotkreuz, Switzerland

³Institute for Pharmacology and Preventive Medicine, Menzelstrasse 21, 15831 Mahlow, Germany

Corresponding author: Claudia Zemmrich, claudia.zemmrich@ippmed.de

Abstract

Background:

Immunochromatographic rapid antigen tests (RATs) emerged onto the COVID-19 pandemic testing landscape to aid in the rapid diagnosis of people with suspected SARS-CoV-2 infection. RATs are particularly useful where RT-PCR is not immediately available and symptoms suggestive of a high viral load and infectiousness are assumed. Several lateral flow immunoassays have been authorized for use under EUA and/or the CE mark, presenting varying overall clinical performance data generated by the manufacturer or by independent investigators. To compare the real-world clinical performance of commercially available rapid chromatographic immunoassays intended for the qualitative detection of SARS-CoV-2, we performed a systematic meta-analysis of published data.

Methods:

We searched the MEDLINE®, Embase, BIOSIS and Derwent Drug File (ProQuest) for manufacturer-independent prospective clinical performance studies comparing SARS-CoV-2 RATs and RT-PCR assays. Only studies on lateral flow assays not needing a separate reader for retrieving the result were included, if data were available on viral load, patients' symptom status, sample type, and PCR assay used. For better data comparability, recalculation of the studies' single performance data confidence intervals using the exact Clopper–Pearson method was applied.

Results:

We could include 19 studies (ten peer-reviewed) presenting detailed clinical performance data on 11,209 samples with 2449 RT-PCR-positives out of study prevalence rates between 1.9–100 % and between 50–100% symptomatic samples. Four studies directly compared two to three different RATs and 15 studies compared one RAT to RT-PCR. Overall specificity ranged, with one test outlier, between 92.4% (87.4–

NOTE: This preprint reports new research that has not been certified by peer review and should not be used to guide clinical practice.

95.9) and 100% (99.7–100), and overall clinical sensitivity varied between 28.9% (16.4–44.3) and 98.3% (91.1–99.7), depending on assay, population characteristics, viral load, and symptom status. Sensitivity in high-viral-load samples (cycle threshold ≤ 25) showed a considerable heterogeneity among the assays ranging from 66.7% to 100%.

Conclusion:

Only two RATs offered sufficient manufacturer-independent, real-world performance data supporting use for the detection of current SARS-CoV-2 infection in symptomatic or high-viral-load patient populations. Reliable positive predictive values require testing of symptomatic patients or asymptomatic individuals only in case of a high pre-test probability. If RATs are used for screening of asymptomatic cases in low-prevalence scenarios, a lower positive predictive value of the result has to be considered.

Background

Real-time reverse transcriptase PCR (RT-PCR) is the current gold standard to detect acute SARS-CoV-2 infection using nasopharyngeal, nasal, or oropharyngeal swab samples containing a measurable amount of viral RNA. It requires a professionally run laboratory with molecular-biological competence, as well as transport infrastructure between the place of sample collection and the lab. Immunochromatographic rapid antigen tests (RATs) emerged on the pandemic testing landscape to support the rapid diagnosis of individuals suspected of SARS-CoV-2 infection, either by presenting symptoms or by contact with cases. These point-of-care tests are less clinically sensitive than RT-PCR tests but overall offer a reliable specificity. Several lateral flow immunoassays have been authorized for use under EUA and/or the CE mark,^{1,2} presenting clinical performance data generated mostly by the manufacturer and containing heterogeneous patient populations. The Foundation for Innovative New Diagnostics (FIND) collates new diagnostic SARS-CoV-2 test developments and manufacturer-independent validation studies of newly commercialized tests on a regularly updated webpage.³

Sensitivity and specificity of RATs depend on numerous variables contributing to the final test results. Their lower limit of detection is at a higher viral load compared with RT-PCR tests.⁴ Most RATs are intended for use in patients up to 5 to 7 days after symptom onset, to increase the probability of having sufficiently high viral load for detection. If RATs are used to assess asymptomatic contact samples, time from symptom onset is not available and can only be assumed by the date of the index case contact. The use of a RAT for screening within a low-prevalence population may not be appropriate, as fewer cases with a detectable high viral load are expected within this group, decreasing the positive predictive value of the test accordingly.⁵

Direct assay comparison studies minimize further heterogeneity parameters (e.g. PCR assay performance differences or varying circulating virus cell characteristics) but are limited due to the invasive character of repeat sample extraction and a high number of required screened persons to detect a sufficient number of positive cases.

In order to provide clarity on the real-world clinical performance of RATs, we compiled all available manufacturer-independent, prospectively collected, clinical data using rapid chromatographic immunoassays, which were commercially available as of November 20, 2020 and intended for the qualitative detection of SARS-CoV-2 present in human nasopharynx, in individuals suspected of SARS-CoV-2 infection. We aimed to harmonize the data regarding the aforementioned performance-impacting factors as much as possible using mathematical methods, to ensure that the data are to a minimum comparable, despite varying methods of presentation in the publications considered for this analysis.

Materials and Methods

Search strategy

We searched MEDLINE®, EMBASE, BIOSIS, and Derwent Drug File at Host "ProQuest" for any clinical performance studies using a commercial SARS-CoV-2 RAT for the following search terms: "COVID-19" OR "SARS-CoV-2" OR "2019-nCoV" OR "coronavirus disease 2019" OR "novel corona virus" OR MESH Entries for Coronaviridae (incl. narrow terms) OR EMTREE Entries for Coronaviridae (incl. narrow terms) OR MESH/EMTREE Entries for "severe acute respiratory syndrome" (incl. narrow terms) AND "rapid antigen test*" OR "rapid antigen assay*" OR "standard Q covid-19 ag" AND "sensitivity" OR "specificity" OR "clinical performance" OR "positive agreement" OR "negative agreement" OR "concordance" OR "validation" OR "evaluation" OR "accuracy". After exclusion of 40 findings not representing clinical trials (29 findings) or not meeting our inclusion criteria (11 findings), eight clinical trials were extracted, along with two further peer-reviewed trial manuscripts that were found manually. Secondly, we searched the studies listed on the Foundation for Innovative Diagnostics (FIND) website (only final, no interim study reports were considered), the European Commission (EC) COVID-19 Diagnostic Devices and Test Method database website (both accessed December 10, 2020), and the diagnosticsglobalhealth.org website (accessed December 14, 2020); nine pre-print clinical RAT studies were added to our finding list. A summary of the search result is presented in Figure 1.

Selection criteria

Only studies investigating clinical performance data of RATs not requiring a separate reading device compared to RT-PCR were considered, if they were performed independent of funding by the manufacturer or distributor. Retrospective laboratory studies were excluded, information on cycle threshold (Ct) values and symptom status had to be available. To ensure comparability, only nasopharyngeal or combined oro- and nasopharyngeal sample types were considered, either directly performed at the point-of-care or at a laboratory after sample transport in viral transport media. If stated, the target gene(s) of the RT-PCR assay was/were reported. No further exclusion criteria on the studies considered for this analysis have been defined. Parameters such as different classification of a negative sample result according to Ct >38 or >40, or the definition of symptoms qualifying as a "symptomatic patient suspected of COVID-19" varied between the studies.

The data were extracted to an electronic database and stratified according to the testing devices. Clinical performance data per Ct group, and stratification for symptoms versus no symptoms present, were calculated from the single data in case they were available. Finally ten peer-reviewed studies and nine pre-print manuscripts could have been included in the meta-analysis, investigating five different RATs:

the SARS-CoV-2 Rapid Antigen Test (Roche), equaling the STANDARD Q COVID-19 Ag Test (SD Biosensor, henceforth called “Roche/SDB”); the Panbio™ COVID-19 Ag Test (Abbott, henceforth called “Abbott”); the COVID-19 Ag Respi-Strip® (Coris Bioconcept, henceforth called “Coris”); the COVID-Viro® (AAZ-LMB, henceforth called “AAZ”); and the BIOCREDIT COVID-19 Ag (RapiGEN, henceforth called “RapiGEN”).

Statistical analyses

Performance results from the comparison of the RATs against the RT-PCR methods are reported as sensitivity and specificity. If available, sensitivity is also provided for cohorts dependent on the Ct values. As confidence intervals reported in the manuscripts used different methods, all confidence intervals have been recalculated for our review using the exact Clopper–Pearson method to allow for better comparability.

The meta-analysis of the performance results of the RATs against the RT-PCR methods was performed for the Roche/SDB, Abbott, and Coris RATs using the statistical software R.

The metaprop function from the “meta” package was used to calculate the effect size for each different test and overall. The results of the AAZ and RapiGEN RATs are also included, even though only one study was available for those tests. The results are shown as a forest plot summarizing the sensitivities found in the different studies.

The bivariate model of Reitsma et al.⁶ was fitted as a linear mixed model and variance components are estimated by restricted maximum likelihood (REML), using the reitsma function from the “mada” package for each system investigated in more than one study. The results are presented as a summary receiver operating characteristic (SROC) curve plot⁷ including the results of all systems (also those investigated in only one study). The single studies, summary estimates, SROC curves, and confidence regions are depicted.

The relationship between sensitivity and viral load represented by Ct value is visualized in a scatterplot. The single study results for sensitivity below a certain Ct-threshold are plotted against these Ct values categorized by the different systems. If in a single study sensitivity estimates for more than one Ct value were available, those are connected by a line.

Results

Study population characteristics

Altogether, we included data from 19 clinical studies providing single data on 11,209 samples, including 2449 samples with confirmed SARS-Cov-2 by RT-PCR, see Tables 1–4. Three studies directly compared

the Roche/SDB RAT with the Abbott,⁸ Coris,⁹ or RapiGEN¹⁰ RATs, respectively, and one study compared the Roche/SDB RAT with both the Abbott and the AAZ RAT,¹¹ all against RT-PCR. 15 studies compared the clinical performance of one RAT against RT-PCR: six studies compared the Roche/SDB;¹²⁻¹⁷ six compared the Abbott;¹⁸⁻²³ and three compared the Coris RAT.²⁴⁻²⁶ All studies provided descriptions of the study populations, regarding mean age and gender distribution (data not shown), as well as symptoms, prevalence rates, and Ct of the RT-PCR-positive samples. Prevalence rate – here meaning the number of RT-PCR-positive samples within the overall study population – varied between 1.9–100%. The prevalence of SARS-CoV-2 in some of these studies did not reflect the prevalence in the local populations, as additional pre-specified testing criteria qualified patients for study entry and created a preselection bias.

Sample characterization according to Ct value of the RT-PCR-positive samples showed a mean Ct value between 20–33. The definition of high viral load varied considerably within the trials, from Ct <18.577 in an Egyptian study¹⁰ and <37 in a Ugandan trial.¹⁷ The Ct groups, which were summarized and analyzed by group, also differed considerably. A majority, but not all, have chosen to report groups of Ct ≤20, ≤25, ≤30, and ≤35. The cobas® (Roche) and Allplex / Seegene® were the most frequently used RT-PCR assays, but even within a single study up to five different assays were used.⁹ Mostly the E-gene was detected, sometimes the RdRp-gene or N-gene as well, and sometimes the target was not stated. One study did not report Ct values but reported RNA copies/mL.¹¹

Local definitions of “patients suspected of SARS-CoV-2 infection” either included only patients presenting clinical symptoms, or also asymptomatic persons with recent direct contact with suspected or confirmed cases. Some asymptomatic contact case groups were limited to healthcare professionals with recent contact, but others included travelers returning from risk areas, which meant that this population varied in terms of pretest probability. The proportion of symptomatic versus asymptomatic samples highly varied, ranging between 14.0–97.8%. Not every report clearly stated the ratio between symptomatic RT-PCR-positive and symptomatic RT-PCR-negative samples. This ratio seemed to differ considerably; some studies reported overall >90% symptomatic persons gaining approximately 15% RT-PCR-positive samples, and the study which investigated only 14% symptomatic persons tested 34.4% RT-PCR-positive samples.

Clinical performance of the RATs

Overall specificity of the investigated RATs ranged between 92.4–100% with one outlier (45%),¹⁰ see Tables 1–4. The two rapid test assays with the most comprehensive available database of >6 studies, the Roche/SDB and Abbott, reported a specificity of >97% in the majority of the according trials, meeting the requirement threshold defined by health authorities.^{5,27} The Coris RAT ranged between 95.8–100% for

specificity but this was combined with unacceptable low sensitivity. The AAZ showed very good results with specificity of 100% and sensitivity of 84.1% but was only evaluated in a single study and RapiGEN showed unacceptable low specificity of 45% in the only published study, combined with low sensitivity as well.

Pooled reported overall sensitivities ranged between 37.7–84.1% (95% CI 27.8–48.7 and 77.0–89.3 respectively), demonstrating considerable heterogeneity between the assays but better comparability within the assays, see Figure 2. Only the AAZ, Roche/SDB and Abbott tests showed a comparable overall performance and comply with Health Authorities' requirement of a sensitivity >80%,^{5,27} see according SROC plots in Figure 3. Expectedly, all assays performed better in samples with high viral loads, but sensitivity dropped more rapidly at Ct >20 for the Coris and RapiGEN tests and less rapidly for the Roche/SDB and Abbott tests, see Figure 4. All assays showed a low sensitivity at Ct >30–32 and variable accuracy at Ct between 25–30.

Critical appraisal of data comparability

Our review presents an overview of available manufacturer-independent, clinical performance data of commercial SARS-CoV-2 RATs not requiring a reading instrument. Altogether, 19 studies investigating five different RATs presented detailed population characteristics and Ct values. A rather comprehensive database exists for two assays, while the other three tests present insufficient data and do not meet the WHO requirement thresholds for sensitivity and specificity.

There is one major concern to be highlighted when comparing performance data of RATs originating from different performance studies:

Data presented in different trials must not directly be compared with each other as numerous variables impact on the resulting performance values. A direct comparison can only be made within one trial, performing a head-to-head comparison of different tests within an identical setting or when using several studies such that effects average out. Several variables impact considerably on clinical performance numbers and prohibit any direct comparison:

1. Preanalytical influencers
 - a. Sample type and way of sampling (same or opposite nostril, combined oro-/nasopharyngeal vs nasopharyngeal sample only, order of sampling)
 - b. Collection device, transport media and volume vs direct testing without dilution by transport media
 - c. Time to test and storage/transport conditions, time delay before processing

2. Analytical influencers

- a. Viral load of the sample and viral load distribution in the respective cohort, represented by Ct or RNA copies/mL
- b. Analytical sensitivity and specificity of the assay
- c. PCR assay specifics, different target genes (E-/RdRp-/N-gene etc.)
- d. Across-laboratory differences (e.g. definition of a positive sample starting at Ct <38 or <40)

3. Clinical parameters of the tested subject

- a. Days post symptom onset of sampling
- b. Asymptomatic vs symptomatic status, definition of symptoms “suspect of SARS-Cov-2 infection”
- c. Severity of symptom presentation

It must be noted that symptom classifications considerably differed between the studies. One study even investigated different populations for self-reported versus physician-defined symptoms. A uniform definition of “clinical symptoms suggestive for SARS-CoV-2 infection” would be desirable but is currently not available.

Due to methodological reasons, the detection limit for SARS-CoV-2 RNA material out of clinical samples tested by RT-PCR is always lower than the detection limit for SARS-Cov-2 antigen. With decreasing viral load, the detectability of even the best performing RAT deteriorates. Cell culture studies on discrepant test results (positive on RT-PCR, negative on RAT) show rare to no virus growth at such samples with high Ct values.^{9,28} This translates into very limited to no infectiousness of the infected patients, even if RT-PCR may still show positive signals for up to three more weeks after peak Ct value.^{29,30}

Conclusion

Only two RATs offer sufficient real-world and manufacturer-independent performance data to recommend their use for the detection of current SARS-CoV-2 infection in symptomatic or suspected high-viral-load patient populations. A high positive predictive value requires testing at a high pre-test probability setting. This requires careful preselection and confirmation of recent contact to confirmed cases and/or knowledge about the underlying local population prevalence. If RATs are used for screening of asymptomatic cases in low-prevalence scenarios, a lower positive predictive value of the according result has to be considered.

Acknowledgements

Editorial assistance was provided by Elements Communications Ltd, Westerham, UK and funded by Roche Diagnostics.

Competing interests

All authors have completed the ICMJE uniform disclosure form at www.icmje.org/coi_disclosure.pdf and declare support from Roche Diagnostics for the submitted work. Johannes Hayer and Dusanka Kasapic are employees of Roche Diagnostics. Claudia Zemmrich works as a freelance contractor for Roche Diagnostics.

Data availability

The data supporting this meta-analysis are from previously reported studies and datasets which have been cited. The processed data are available from the corresponding author upon request.

Funding

This work was supported by Roche Diagnostics.

Figures

Figure 1. Flowchart of search results

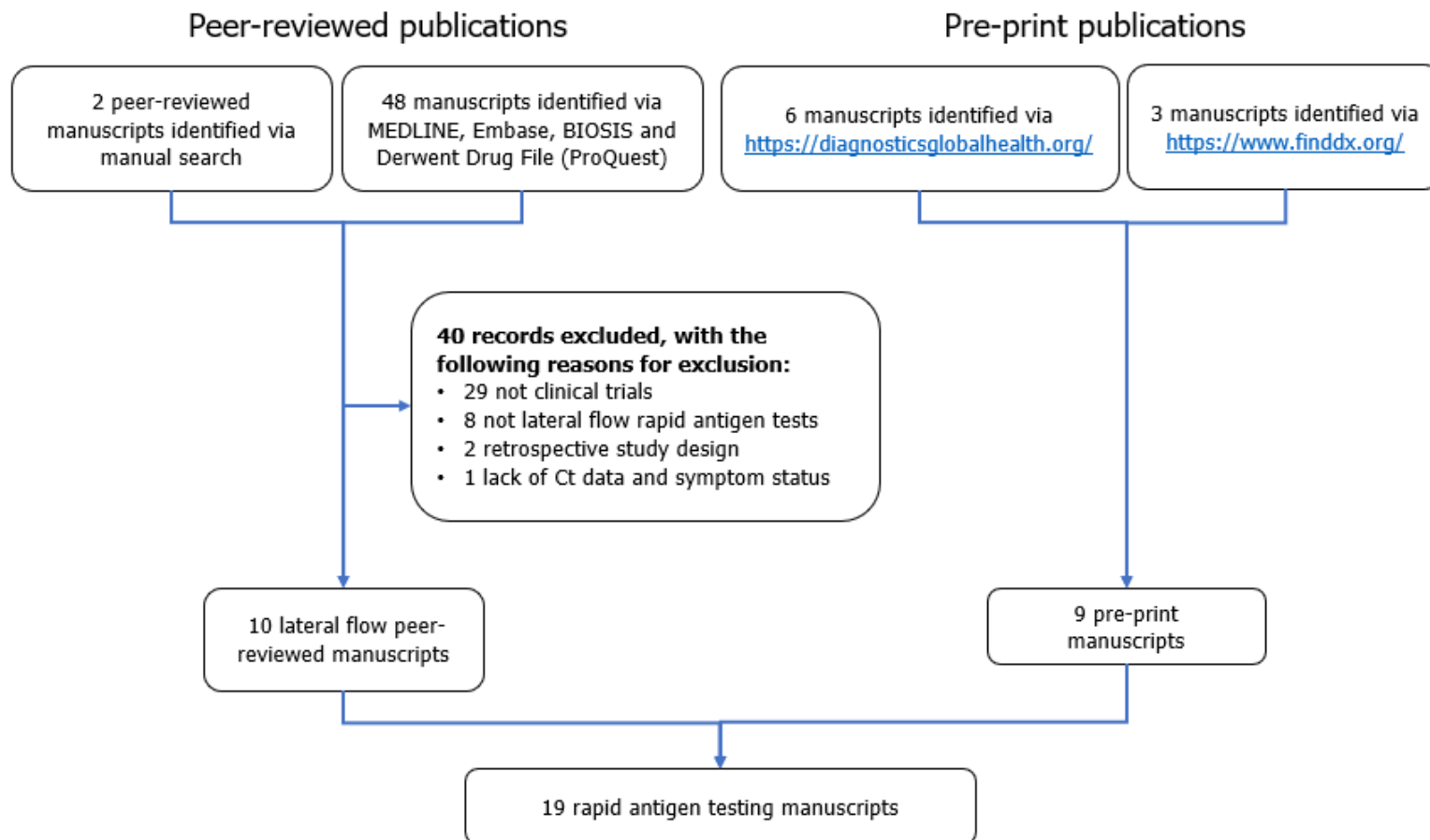


Figure 2. Forest plot of studies evaluating rapid antigen test sensitivity, grouped by test



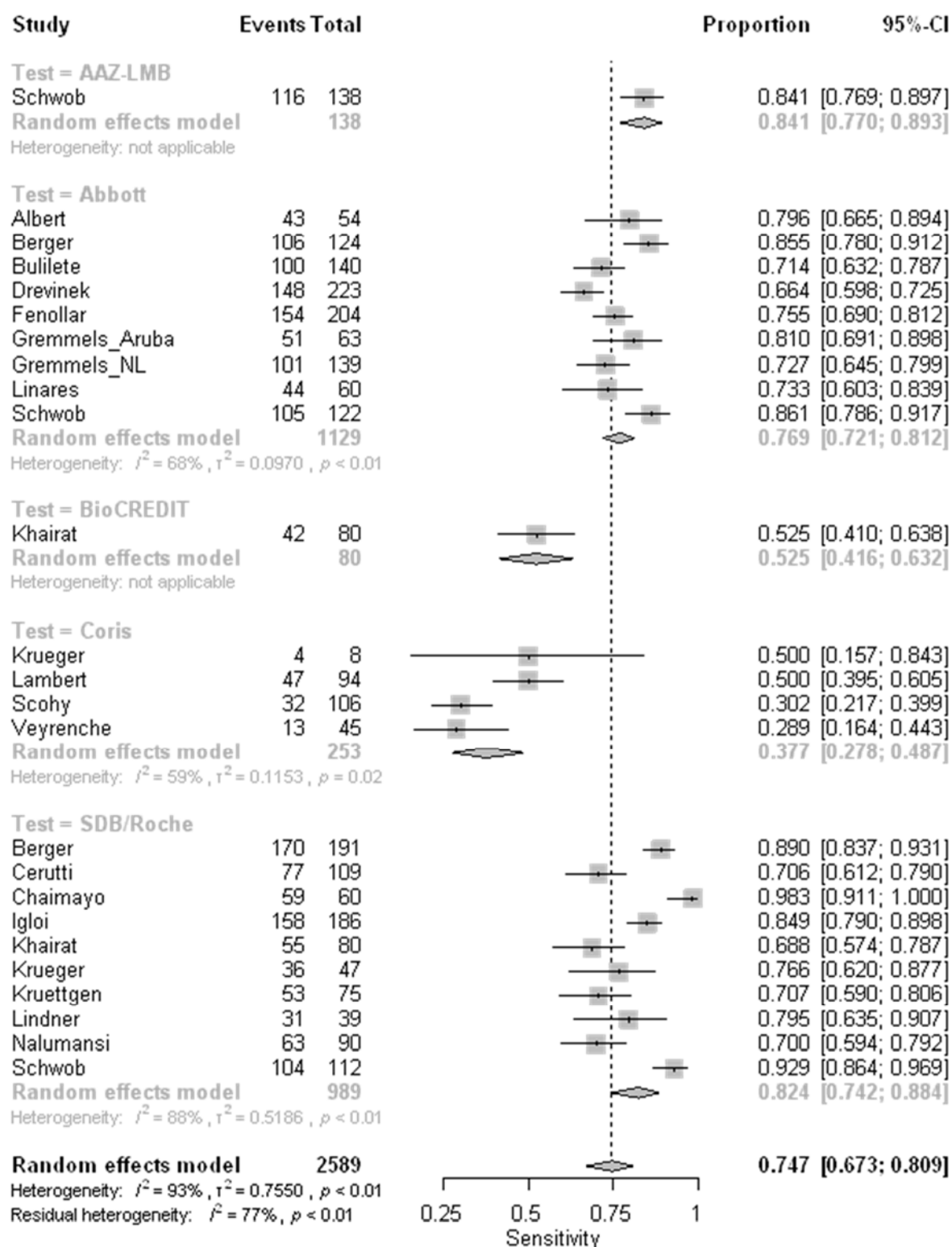


Figure 3. SROC plots for all*

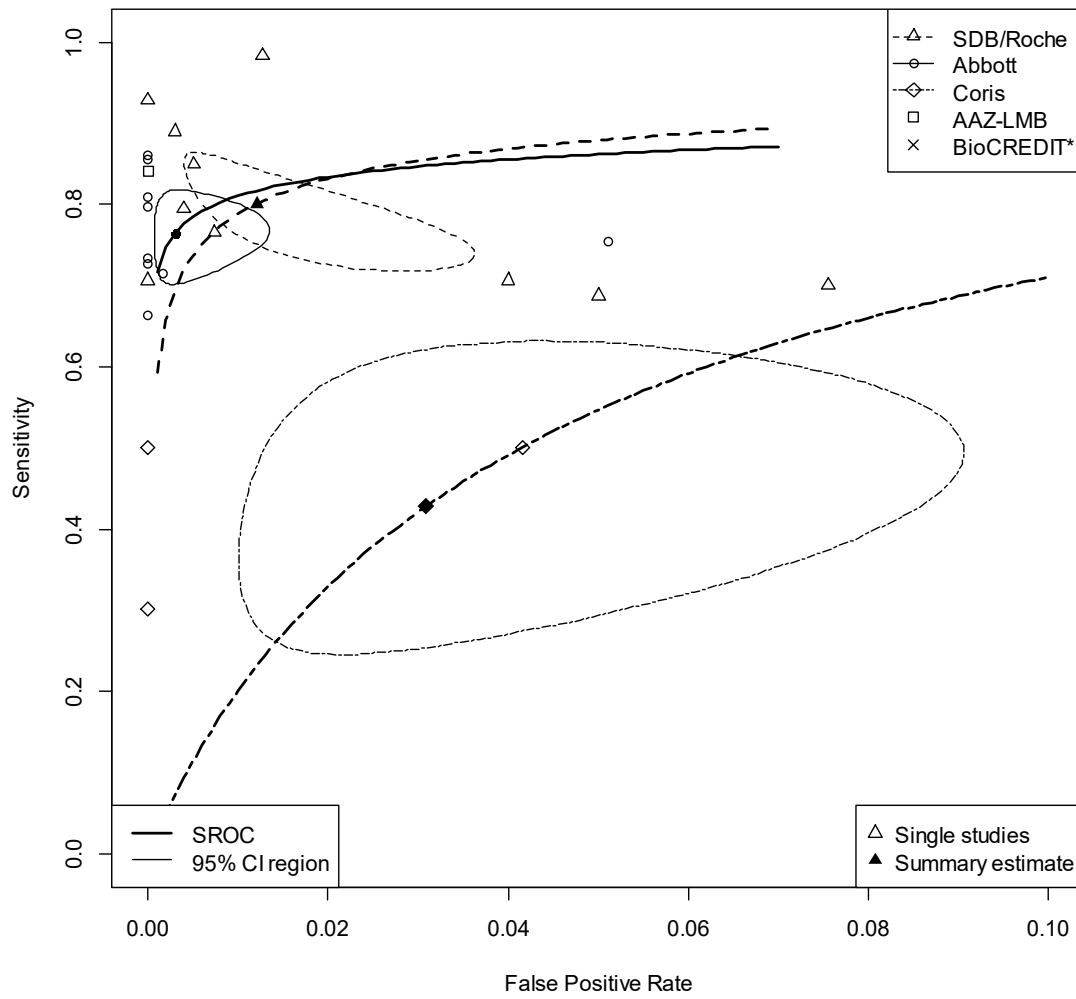

*Study result for BioCREDIT lies outside the plotting region ($x=0.55$, $y=0.525$) and is therefore not shown.

Figure 4. Rapid antigen test sensitivity according to viral load, represented by Ct value groups

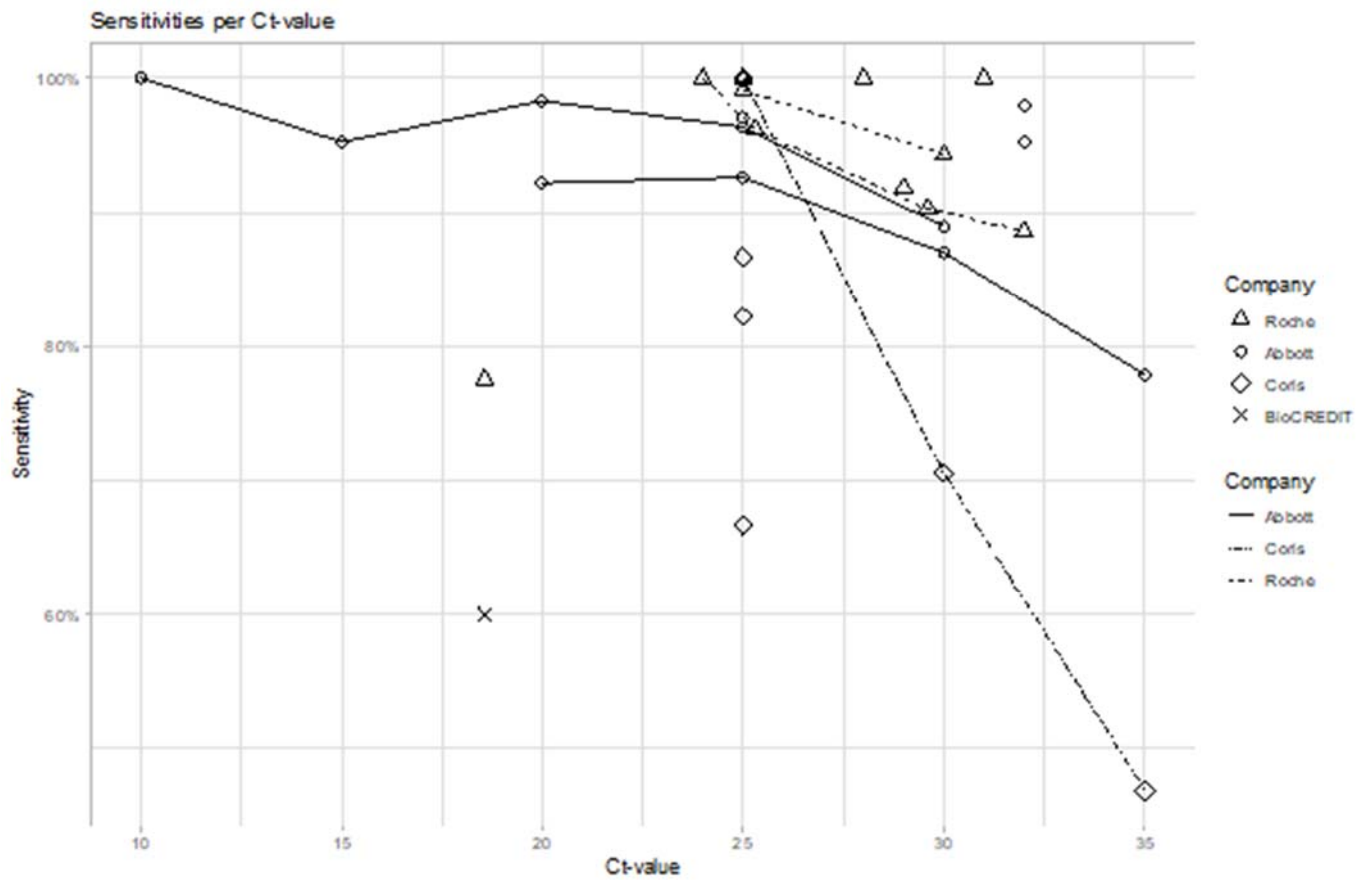

Tables

Table 1. Study population characteristics and overall performance: Roche/SDB only

	Chaimayo¹²	Igloi¹⁵	Lindner¹⁶	Cerutti¹³	Krüttgen¹⁴	Nalumansi¹⁷
Country	Thailand	Netherlands	Germany	Italy	Germany	Uganda
Participants (N)	454	970	289	330	150	262
PCR+ (n)	60	186	39	109	75	90
Prevalence (%)	13.2	19.2	13.5	33	50	34.4
Symptomatic (%)	95.0 (PCR+)	91.3 (overall)	97.6 (overall)	56.1 (overall); 95.4 (PCR+)	Not reported	14 (PCR+)
PCR manufacturer and assay	Seegene Allplex™ 2019-nCoV Assay	Roche cobas® SARS-CoV-2 Assay	Roche cobas® SARS-CoV-2 Assay; TIB MOLBIOL SARS-CoV-2 E-gene®	Roche cobas® SARS-CoV-2 Assay; Seegene Allplex™ 2019-nCoV Assay; DiaSorin Simplexa®	Altona Diagnostics RealStar® SARS-CoV-2	Metabion custom oligo nucleotides
Specimen Ct values (min – max); Mean/Median	10.5–39.0; Mean 22.8±6.7 (E-gene)	15.6–37.4; Median 23.6	17.3–35.5; Mean 23.7±5.5; Median 21.9	12.3–38.1	<20–≥ 35	≤29–39
Antigen test performance						
Overall specificity (95% CI)[#]	98.7 (97.1–99.5)	99.5 (98.7–99.9)	99.6 (97.8–100)	100 (98.3–100)	96.0 (88.8–99.2)	92.4 (87.4–95.9)
Overall sensitivity (95% CI)[#]	98.3 (91.1–99.7)	84.9 (79.0–89.8)	79.5 (63.5–90.7)	70.6 (61.2–79.0)	70.7 (59.0–80.6)	70 (59.4–79.2)
Sensitivity by Ct (95% CI)[#]	Ct≤31: 100	Ct≤25: 99.1 (95.2–100); Ct≤30: 94.3 (89.6–97.0)	Ct≤24: 100 (85.2–100); Ct≤25.3: 96.2 (80.4–99.9); Ct≤29.6: 90.3 (74.2–98.0); Ct≤32: 88.6 (73.3–96.8)	Ct<28: 100; Ct 28–30: 38.5; Ct 30–35: 26.7; Ct>35: 9.1	Ct<25:100; Ct 25–<30: 95.0; Ct 30–<35: 44.8; Ct>=35: 22.2	Ct≤29: 91.9% (78.1–98.3); Ct 30–37: 54.5%; Ct 38–39: 55.6%

[#]Values were recalculated using the original data and confidence intervals were calculated using the exact Clopper–Pearson method.

Table 2. Study population characteristics and overall performance – Panbio, Abbott only

	Albert²³	Linares²²	Drevinek¹⁸	Bulilete¹⁹	Fenollar²⁰	Gremmels²¹	
Country	Spain	Spain	Czech Republic	Spain	France	The Netherlands	Aruba
Participants (N)	412	255	591	1369	341	1367	208
PCR+ (n)	54	60	223	140 (7 inconclusive)	204	139	63
Prevalence (%)	13.1	23.5	37.7	10.2	59.8	10.2	30.3
Symptomatic (%)		72.2	49.1 (overall); 75.3 (PCR+)	49.7 (overall); 62.1% (PCR+)	53.4 (overall)	97.3 (overall)	Not reported
PCR manufacturer and assay	ThermoFisher TaqPath™ COVID-19 Combo Kit	Seegene Allplex™ 2019-nCoV Assay	Seegene Allplex™ 2019-nCoV Assay	ThermoFisher TaqPath™ COVID-19 kit	Credo Diagnostics VitaPCR SARS-CoV-2 assay	Seegene Allplex™ 2019-nCoV Assay	Seegene Allplex™ 2019-nCoV Assay
Specimen Ct values (min – max); Mean/Median	<10 to >30	(0–37.8); Median 23.3	<10 to ≥35	Mean 20.3±6.5 (N-gene); Mean 21.9±6.5 (S-gene); Mean 21.0±6.7 (ORF-gene)	<10–34	Mean 27.5±6.0 (N-gene); Mean 24.7±5.7 (E-gene); Mean 26.4±5.6; (RdRp-gene)	Not reported
Antigen test performance							
Overall specificity (95% CI)[#]	100 (99.0–100)	100 (98.1–100)	100 (99.0–100)	99.8 (99.4–100)	94.9 (89.8–97.9)	100 (99.7–100)	100 (97.1–100)
Overall sensitivity (95% CI)[#]	79.6 (66.5–89.4)	73.3 (60.3–83.9)	66.4 (59.8–72.5)	71.4 (63.2–78.7)	75.5 (69.0–81.2)	72.7 (64.5–79.9)	81.0 (69.1–88.8)
Sensitivity by Ct (95% CI)[#]	Ct<25: 100	Ct<25: 97.1 (84.7–99.9); Ct 25–30: 77.8; Ct 30–35: 30.0; Ct 35–40: 14.0	Ct<20: 92.2 (81.1–97.8); Ct<25: 92.6 (86.3–96.5); Ct<30: 87.0 (80.8–91.7); Ct≤35: 77.9 (71.3–83.6)	Not reported	Ct<10: 100; Ct<15: 95.2; Ct<20: 98.3; Ct<25: 96.4; Ct<30: 89.0	Ct<32: 95.3 (89.3–98.5)	Ct<32: 98.0 (89.1–99.9)

[#]Values were recalculated using the original data and confidence intervals were calculated using the exact Clopper–Pearson method.

Table 3. Study population characteristics and overall performance – CORIS RESPI-strip only

	Lambert-Niclot ²⁴	Veyrenche ²⁵	Scohy ²⁶
Country	France	France	Belgium
Participants (N)	138	45	148
PCR+ (n)	94	45	106
Prevalence (%)	68.1	NA (only positive clinical samples)	71.6
Symptomatic (%)	Not reported	100	88.5 (overall)
PCR manufacturer and assay	Altona Diagnostics RealStar® SARS-CoV-2; Anatolia Geneworks Bosphore novel coronavirus (2019 nCoV) detection kit; Roche cobas® SARS-CoV-2 Assay; Seegene Allplex™ 2019-nCoV Assay	Seegene Allplex™ 2019-nCoV Assay	Primerdesign genesig® Real-Time PCR coronavirus COVID-19
Specimen Ct values (min – max); Mean/Median	≤10 to >40	<15 to >40	(16–38); Mean 31.4; Median 33
Antigen test performance			
Overall specificity (95% CI)[#]	100 (92.0–100)	NA	100 (91.6–100)
Overall sensitivity (95% CI)[#]	50.0 (39.5–60.5)	28.9 (16.4–44.3)	30.2 (21.7–39.9)
Sensitivity by Ct (95% CI)[#]	Ct<25: 82.2	Ct≤25: 86.7 (59.5–98.3); Ct>25: 0.0 (0.0–11.6)	Ct<25: 100; Ct<30: 70.6; Ct<35: 46.9

[#]Values were recalculated using the original data and confidence intervals were calculated using the exact Clopper–Pearson method.

Table 4. Study population characteristics and overall performance – studies assessing multiple antigen assays

	Schwob ¹¹			Berger ⁸		Krüger ⁹		Khairat ¹⁰	
Country	Switzerland			Switzerland		Germany + UK		Egypt	
Antigen assay	Roche/SDB	Panbio, Abbott	COVID-Viro AAZ-LMB	Panbio, Abbott	Roche/SDB	Roche/SDB	CORIS RESPI-strip	Roche/SDB	BioCREDIT
Participants (N)	928			1064		1688*		100	
Participants per assay (n)	333	271	324	535	529	1263	425	100	100
PCR+ (n)	112	122	138	124	191	47	8	80	80
Prevalence (%)	33.6	45.0	42.6	23.2	36.1	3.7	1.9	75	75
Symptomatic (%)	96 (overall)			97.8 (PCR+)		84.4 (overall)	68.9 (overall)	Unreported	
PCR manufacturer and assay	Roche cobas® SARS-CoV-2 Assay; in-house RT-PCR			Roche cobas® SARS-CoV-2 Assay		Germany: TIB MOLBIOL SARS-CoV-2 E-gene; Abbott RealTime SARS-COV-2 assay; Roche cobas® SARS-CoV-2 Assay; Seegene Allplex™ 2019-nCoV Assay UK: Primerdesign genesig® Real-Time PCR coronavirus COVID-19		Unreported	
Specimen Ct values (min – max); Mean/Median	NA (Results per viral load available)			(14.2–39.7); Mean 22.5±5.1; Median 21.5		NA	NA	Median 18.57	
				(14.2–39.7); Mean 22.4±5.4; Median 21.8	(14.4–37.4); Mean 22.6±4.9; Median 21.0				
Antigen test performance									
Overall specificity (95% CI) [#]	100 (98.3–100)	100 (97.6–100)	100 (98.0–100)	100 (99.1–100)	99.7 (98.3–99.9)	99.3 (98.6–99.7)	95.8 (93.4–97.6)	95 (75.1–99.9)	45 (23.1–68.5)
Overall sensitivity (95% CI) [#]	92.9 (86.4–96.9)	86.1 (78.6–91.7)	84.1 (76.9–89.7)	85.5 (78.0–91.2)	89.0 (83.7–93.1)	76.6 (62.0–87.7)	50.0 (21.5–78.5)	68.8 (57.4–78.7)	52.5 (41.0–63.8)

	Schwob ¹¹			Berger ⁸		Krüger ⁹		Khairat ¹⁰	
Sensitivity by Ct (95% CI)[#]						Ct<25: 100 (81.5–100); Ct≥25: 62.1 (42.3–79.3)	Ct<25:66.7 (9.4–99.2); Ct≥25:40.0 (5.3–85.3)	Ct<18.57: 77.5 (61.5–89.2); Ct>18.57: 60.0 (43.3–75.1)	Ct<18.57: 60.0 (43.3–75.1); Ct>18.57: 45.0 (29.3–61.5)

*Total sample number of study was 2417, but one assay was excluded from our analysis for not fulfilling inclusion criteria; [#]Values were recalculated using the original data and confidence intervals were calculated using the exact Clopper–Pearson method;

References

1. US Food and Drug Administration. In vitro diagnostics EUAs. 2020. <https://www.fda.gov/medical-devices/coronavirus-disease-2019-covid-19-emergency-use-authorizations-medical-devices/vitro-diagnostics-euas#individual-antigen>.
2. World Health Organization. Regulation and prequalification. 2020. <https://www.who.int/teams/regulation-prequalification/eul/in-vitro-emergency-use-listing-procedure>.
3. The Foundation for Innovative New Diagnostics. SARS-CoV-2 diagnostic pipeline. 2020. <https://www.finddx.org/covid-19/pipeline/>.
4. Dinnes J, Deeks JJ, Adriano A, et al. Rapid, point-of-care antigen and molecular-based tests for diagnosis of SARS-CoV-2 infection. *Cochrane Database Syst Rev* 2020; **8**: Cd013705.
5. Centers for Disease Control and Prevention. Interim guidance for antigen testing for SARS-CoV-2. 2020. <https://www.cdc.gov/coronavirus/2019-ncov/lab/resources/antigen-tests-guidelines.html>.
6. Reitsma JB, Glas AS, Rutjes AW, Scholten RJ, Bossuyt PM, Zwinderman AH. Bivariate analysis of sensitivity and specificity produces informative summary measures in diagnostic reviews. *J Clin Epidemiol* 2005; **58**(10): 982–90.
7. Rutter CM, Gatsonis CA. A hierarchical regression approach to meta-analysis of diagnostic test accuracy evaluations. *Stat Med* 2001; **20**(19): 286–584.
8. Berger A, Ngo Nsoga MT, Perez-Rodriguez FJ, et al. Diagnostic accuracy of two commercial SARS-CoV-2 Antigen-detecting rapid tests at the point of care in community-based testing centers. *medRxiv* 2020: 2020.11.20.20235341.
9. Krüger LJ, Gaeddert M, Köppel L, et al. Evaluation of the accuracy, ease of use and limit of detection of novel, rapid, antigen-detecting point-of-care diagnostics for SARS-CoV-2. *medRxiv* 2020: 2020.10.01.20203836.
10. Khairat SM, Guindy NEL, Abdel Motaleb MSE, Soliman NS. Evaluation of two rapid antigen tests for detection of SARS-CoV-2 virus. *International Journal of Microbiology and Biotechnology* 2020; **5**(3): 131–4.
11. Schwob JM, Miauton A, Petrovic D, et al. Antigen rapid tests, nasopharyngeal PCR and saliva PCR to detect SARS-CoV-2: A prospective comparative clinical trial. *medRxiv* 2020: 2020.11.23.20237057.
12. Chaimayo C, Kaewnaphan B, Tanlieng N, et al. Rapid SARS-CoV-2 antigen detection assay in comparison with real-time RT-PCR assay for laboratory diagnosis of COVID-19 in Thailand. *Virology* 2020; **17**(1): 177.
13. Cerutti F, Burdino E, Milia MG, et al. Urgent need of rapid tests for SARS CoV-2 antigen detection: Evaluation of the SD-Biosensor antigen test for SARS-CoV-2. *J Clin Virol* 2020; **132**: 104654.
14. Krüttgen A, Cornelissen CG, Dreher M, Hornef MW, Imöhl M, Kleines M. Comparison of the SARS-CoV-2 rapid antigen test to the real star SARS-CoV-2 RT PCR kit. *J Virol Methods* 2020; **288**: 114024.
15. Iglói Z, Velzing J, van Beek J, et al. Clinical evaluation of the Roche/SD Biosensor rapid antigen test with symptomatic, non-hospitalized patients in a municipal health service drive-through testing site. *medRxiv* 2020: 2020.11.18.20234104.
16. Lindner AK, Nikolai O, Kausch F, et al. Head-to-head comparison of SARS-CoV-2 antigen-detecting rapid test with self-collected anterior nasal swab versus professional-collected nasopharyngeal swab. *Eur Respir J* 2020.
17. Nalumansi A, Lutalo T, Kayiwa J, et al. Field evaluation of the performance of a SARS-CoV-2 antigen rapid diagnostic test in Uganda using nasopharyngeal samples. *Int J Infect Dis* 2020.
18. Drevinek P, Hurych J, Kepka Z, et al. The sensitivity of SARS-CoV-2 antigen tests in the view of large-scale testing. *medRxiv* 2020: 2020.11.23.20237198.
19. Bulilete O, Lorente P, Leiva A, et al. Evaluation of the Panbio™ rapid antigen test for SARS-CoV-2 in primary health care centers and test sites. *medRxiv* 2020: 2020.11.13.20231316.
20. Fenollar F, Bouam A, Ballouche M, et al. Evaluation of the Panbio Covid-19 rapid antigen detection test device for the screening of patients with Covid-19. *J Clin Microbiol* 2020.

21. Gremmels H, Winkel BMF, Schuurman R, et al. Real-life validation of the Panbio COVID-19 Antigen Rapid Test (Abbott) in community-dwelling subjects with symptoms of potential SARS-CoV-2 infection. *medRxiv* 2020: 2020.10.16.20214189.
22. Linares M, Pérez Tanoira R, Romanyk J, et al. Panbio antigen rapid test is reliable to diagnose SARS-CoV-2 infection in the first 7 days after the onset of symptoms. *medRxiv* 2020: 2020.09.20.20198192.
23. Albert E, Torres I, Bueno F, et al. Field evaluation of a rapid antigen test (Panbio™ COVID-19 Ag Rapid Test Device) for COVID-19 diagnosis in primary healthcare centres. *Clin Microbiol Infect* 2020.
24. Lambert-Niclot S, Cuffel A, Le Pape S, et al. Evaluation of a rapid diagnostic assay for detection of SARS-CoV-2 antigen in nasopharyngeal swabs. *J Clin Microbiol* 2020; **58**(8).
25. Veyrenche N, Bollore K, Pisoni A, et al. Diagnosis value of SARS-CoV-2 antigen/antibody combined testing using rapid diagnostic tests at hospital admission. *medRxiv* 2020: 2020.09.19.20197855.
26. Scohy A, Anantharajah A, Bodéus M, Kabamba-Mukadi B, Verroken A, Rodriguez-Villalobos H. Low performance of rapid antigen detection test as frontline testing for COVID-19 diagnosis. *J Clin Virol* 2020; **129**: 104455.
27. European Centre for Disease Prevention and Control. Options for the use of rapid antigen tests for COVID-19 in the EU/EEA and the UK. 2020. <https://www.ecdc.europa.eu/en/publications-data/options-use-rapid-antigen-tests-covid-19-eueea-and-uk>.
28. van Beek J, Igloi Z, Boelsums T, et al. From more testing to smart testing: data-guided SARS-CoV-2 testing choices. *medRxiv* 2020: 2020.10.13.20211524.
29. Wölfel R, Corman VM, Guggemos W, et al. Virological assessment of hospitalized patients with COVID-2019. *Nature* 2020; **581**(7809): 465–9.
30. Magleby R, Westblade LF, Trzebucki A, et al. Impact of SARS-CoV-2 viral load on risk of intubation and mortality among hospitalized patients with coronavirus disease 2019. *Clin Infect Dis* 2020.